




AEC-NASA TECH BRIEF

Lawrence Berkeley Laboratory



NASA Tech Briefs announce new technology derived from the U.S. space program. They are issued to encourage commercial application. Tech Briefs are available on a subscription basis from the National Technical Information Service, Springfield, Virginia 22151. Requests for individual copies or questions relating to the Tech Brief program may be directed to the Technology Utilization Office, NASA, Code KT, Washington, D.C. 20546.

Efficient Baffle Prevents Oil Backstreaming in Diffusion Pumps

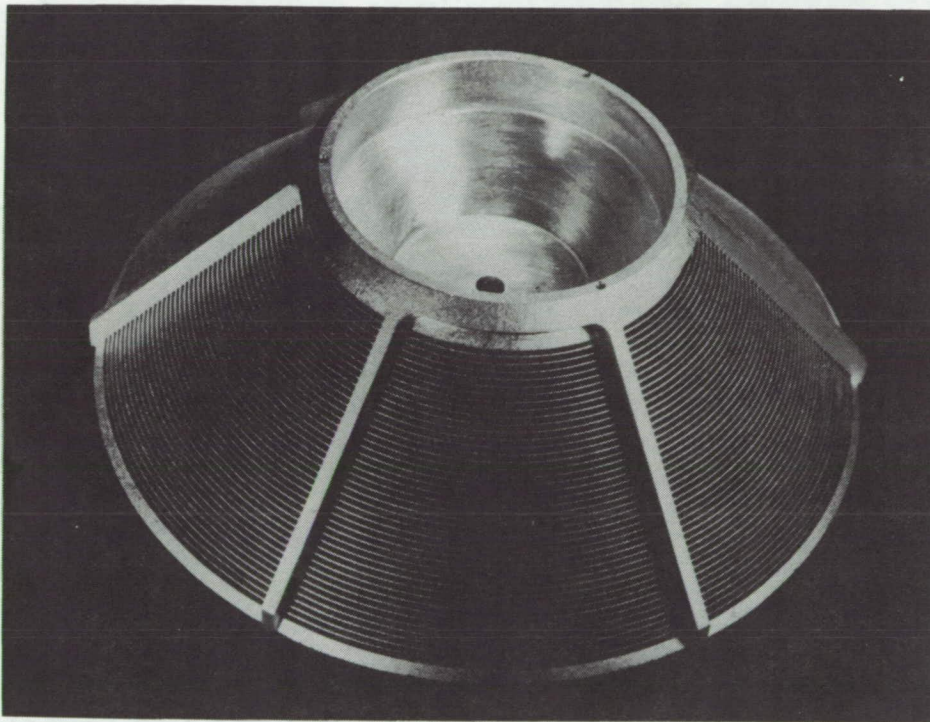

A high-vacuum diffusion pump baffle, positioned immediately above the pump stack, prevents vaporized pumping oil from entering and contaminating the test chamber. The baffle (see figure) consists of a series of concentric rings between the top of the diffusion pump stack and the vacuum pump wall. The pumping vapor is collected and condensed on the baffle and flows back to the pump base along four rods. Although the baffle has no provisions for cooling, test results indicate a

significant reduction in the backstreaming of vaporized pumping oil.

The central cup-shaped section extends downward over the stack, with the baffle section extending upward to the wall of the vacuum pump. Condensed pumping oil, collected on the baffle, flows away from the walls toward the rim of the central cup and then passes down to the reservoir at the pump base.

(continued overleaf)

Note:

Requests for further information may be directed to:
Technology Utilization Division
Lawrence Berkeley Laboratory
University of California
Berkeley, California 94720
Reference: B72-10475

Patent status:

No patent action is contemplated by AEC or NASA.

Source: J. R. Meneghetti
Lawrence Berkeley Laboratory
(LRL-10025)